


# MODUS ESWT RADIAL & MODUS ESWT FOCUSED SHOCKWAVE THERAPY



ESWT SPARE PARTS



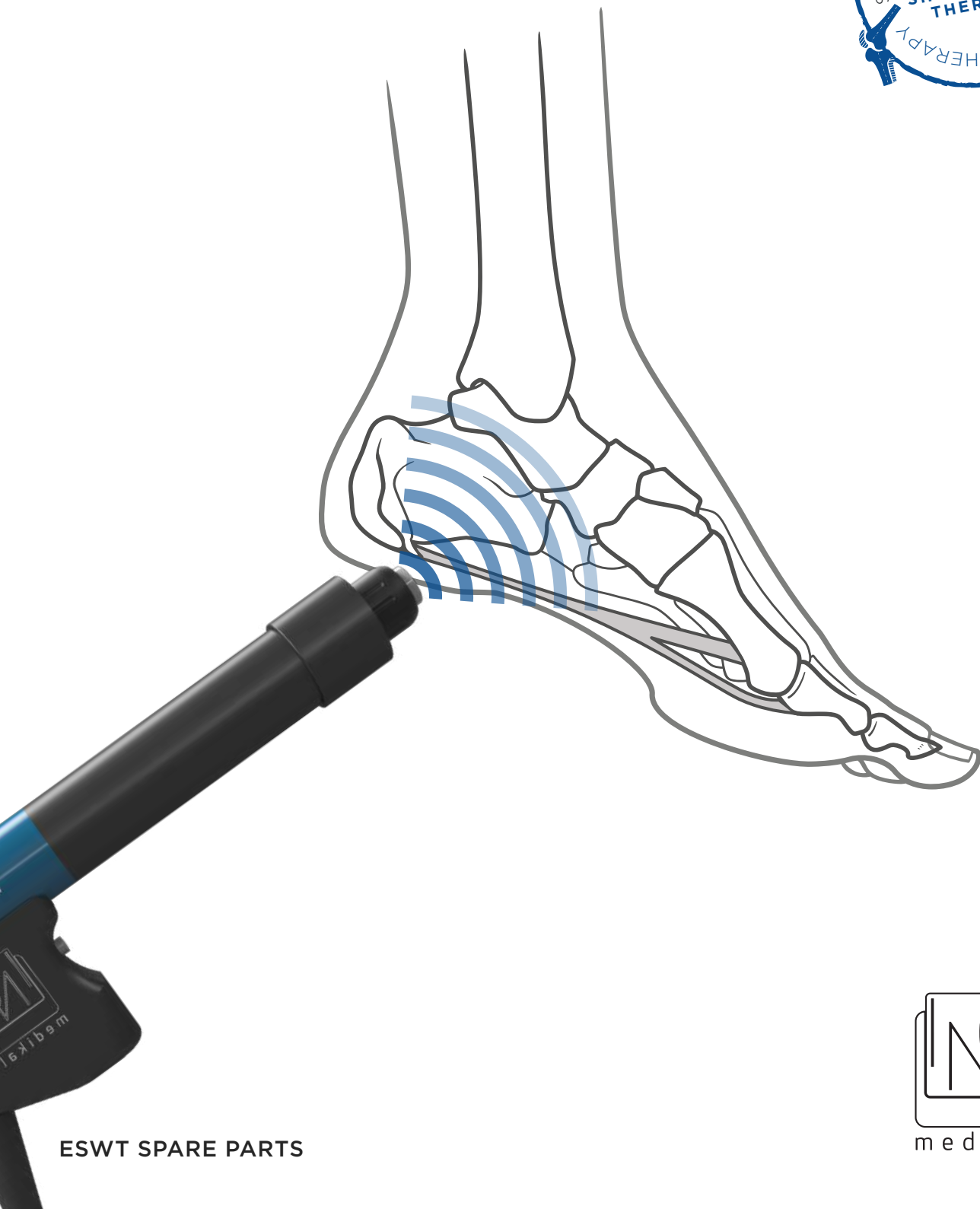

Touch Screen Interface  
With Built-in Compressor  
Modus ESWT Radial

## UNCOVER HEALING POWER MODUS ESWT RADIAL!


Radial shock wave therapy is a multi-purpose treatment system used in orthopedics, physiotherapy, chiropractors, sports medicine, aesthetics, beauty and veterinary medicine.

Radial ESWT devices have a pistol-like handpiece. The bullet in the handpiece moves at a high speed. The bullet hits the tip of the applicator, and the collision creates a pulse of kinetic energy.

In Radial Shock Wave Therapy when shock waves are applied to an anatomical region, they cause mechanical stress and micro-trauma, and this trauma causes the formation of angiogenic factors. Thus, it induces a series of biological reactions that increase blood flow and trigger neovascularization of tissue.

It is an effective treatment method used to eliminate the indications that cause acute and/or chronic pain and to accelerate the recovery without the need for invasive operations-medications.

The shock wave therapy system has no side effects. It does not require medication or surgical intervention in the sessions.

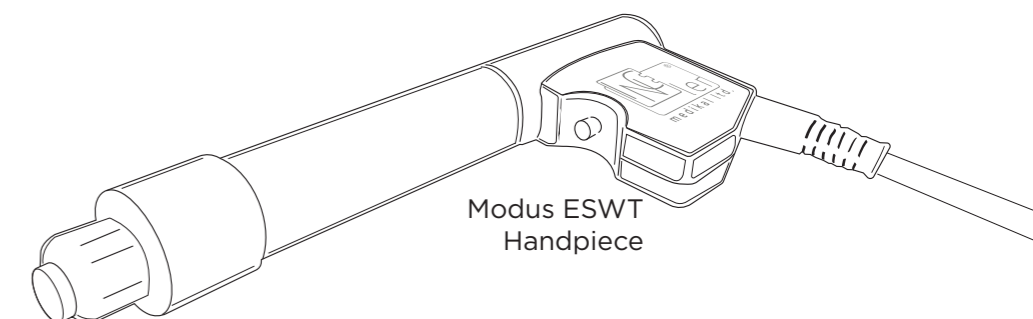

Modus ESWT  
Handpiece

MODUS RADIAL ESWT HANDPIECE  
(WITH ALL APLICATORS)  
Order No.  
IMET-802




## PRODUCT FEATURES

- > As the most powerful Radial Shockwave Therapy device on the market, providing impulses up to 22 Hz and 5 bar
- > Lightweight, portable (10 kg)
- > Opportunity to treat anytime, anywhere
- > The parameters set during the treatment can be easily controlled from the device screen.
- > Ready-made treatment programs and visual, written and video narrative information (Indications and Pathological Applications Guide)
- > Modus ESWT unit with built-in compressor is portable with its original stand
- > It can be easily moved by detaching from its stand when desired.


## INCREASE YOUR APPLICATION EFFICIENCY WITH MODUS HANDPIECES!

- > One-touch on/off feature
- > Suspension system to reduce recoil impact during the treatment
- > Different applicators according to treatment options
- > Easy maintenance feature and easy revision kit replacement
- > Ergonomic design feature that is lightweight and does not tire your hand

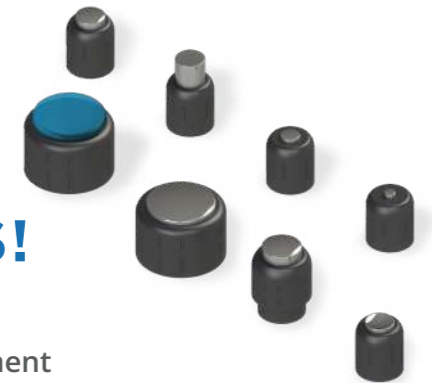

## SOFT APPLICATOR

- > Modus ESWT is the only ESWT device with soft applicator option and this option is making a difference with its unique technology.
- > Soft applicator provides easy usage on painful and sensitive areas.

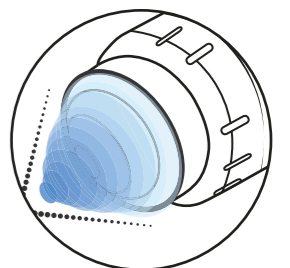



Order No.  
**IMET-889**  
**MODUS RADIAL ESWT HANDPIECE  
(WITH ONLY 15 MM APPLICATORS)**



Order No.  
**IMET-830**  
**MODUS ESWT REVISION KIT**



Order No.  
**IMET-860**  
**EMS COMPATIBLE TYPE 2 REVISION KIT**




Order No.  
**IMET-859**  
**EMS COMPATIBLE TYPE 1 REVISION KIT**




Order No.  
**IMET-861**  
**BTL COMPATIBLE ESWT REVISION KIT**



Order No.  
**IMET-864**  
**STORZ COMPATIBLE ESWT REVISION KIT**

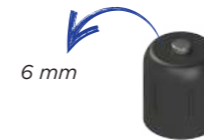

Order No.  
**IMET-805**  
**6 MM APPLICATOR**

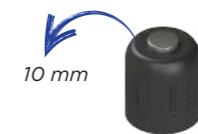

Order No.  
**IMET-806**  
**10 MM APPLICATOR**



Order No.  
**IMET-864**  
**CHATTANOOGA COMPATIBLE ESWT  
ESWT REVISION KIT**



Order No.  
**IMET-864**  
**GYMNA COMPATIBLE ESWT REVISION KIT**




Order No.  
**IMET-807**  
**15 MM APPLICATOR - RADIAL**



Order No.  
**IMET-808**  
**15 MM APPLICATOR - TRIGGER**




Order No.  
**IMET-809**  
**15 MM APPLICATOR - FOCUS**

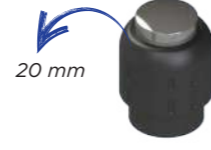

Order No.  
**IMET-810**  
**20 MM APPLICATOR**



Order No.  
**IMET-811**  
**35 MM APPLICATOR**

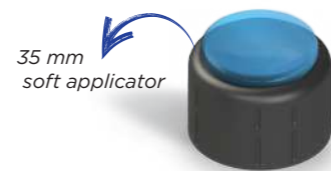

Order No.  
**IMET-812**  
**35 MM MODUS SOFT APPLICATOR**



Order No.  
**IMET-813**  
**35 MM SOFT APPLICATOR COMPATIBLE WITH STORZ**



Order No.  
**IMET-875**  
**POWER CORD**



Order No.  
**IMET-817**  
**MODUS ESWT TROLLEY WITH DRAWER**



Order No.  
**IMET-817-1**  
**MODUS ESWT TROLLEY WITHOUT DRAWER**



Order No.  
**IMPL-728**  
**MODUS ESWT TROLLEY WITH CABINETTE**



Order No.  
**IMET-815**  
**MODUS RADIAL ESWT FILTER (WITH GLASS)**



Order No.  
**IMET-815-1**  
**MODUS RADIAL ESWT FILTER GLASS**



Order No.  
**IMÇ-101**  
**CARRYING BAG**




Touch Screen Interface  
Modus ESWT®  
Focused Shock Waves

## GET ACQUAINTED WITH THE POWER OF MODUS ESWT FOCUSED!

The focusing mechanism directs each shock to a specific point in the body. Separate heads can be used to provide the energy depth. Focused shock waves can target deeper tissues more precisely by propagating all of the energy within the tissue. This is why local pain points, chronic insertional tendinopathies, bone stress injuries, and deep trigger points are best treated with focused shock waves. Patients who cannot tolerate radial treatments for acute injuries can also go to benefit from focused shock wave therapy.

Shock wave therapy is a frequently preferred method in the veterinary medicine, neurology, urology, cardiology, sports medicine, dermatology and aesthetics but especially in the fields of orthopedics and physical therapy.

The non-invasive principle of that the system, does not require anesthesia and surgery. The patient can continue his/her daily routine after the session.

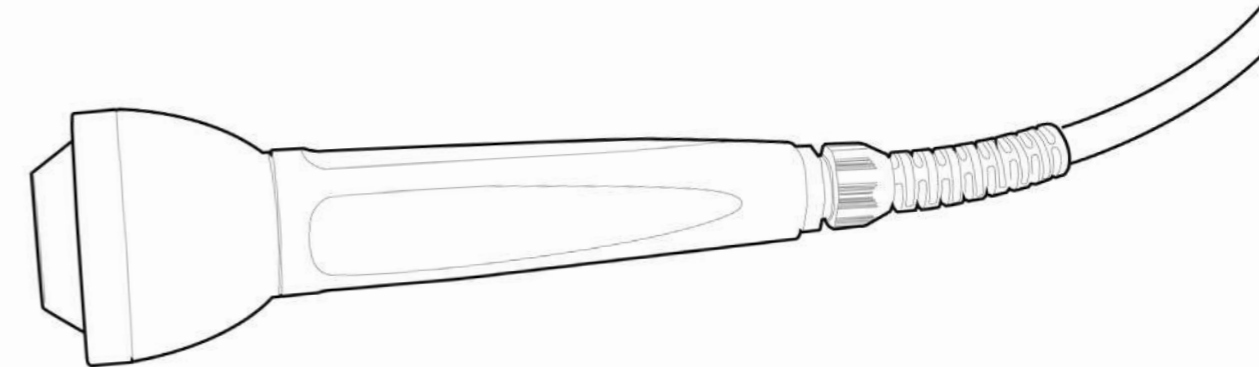

## PRODUCT FEATURES

- > Impulses up to 4 Hz and 0.25mj
- > Head options that provide shock wave transmission with a penetration depth of 0 - 80 mm
- > Power adjustment up to 25 levels
- > Portable, anywhere, anytime treatment
- > Ready-made treatment programs and visual, written and video narrative information(Indications and Pathological Applications Guide)
- > Archive support for the user with patient registration and follow-up menu
- > Ability to change Frequency and Power settings during treatment.
- > The parameters set during the treatment can be easily seen on the device screen and changes can be made during the treatment.
- > Generator and handpiece unit does not require maintenance.
- > Coloured Touchscreen technology provides easy use.

## MODUS FOCUSED HANDPIECE FEATURES

- > In-depth energy delivery to the selected anatomical region
- > Easy and fast use that does not tire the hand;
- > Easy mobility;
- > Maintenance-free handpiece;
- > Head options with easy replacement according to use



Higher Energy Output



High Energy Depth



Fast and Effective  
Treatments



**MODUS FOCUSED ESWT  
HANDPIECE**

Order No.  
IMED-851



**> GIVES MORE  
ENERGY FURTHER**



Order No.

IMED-851

MODUS FOCUSED HANDPIECE  
( WITHOUT HEAD)



Order No.

IMED-864

F-50 HEAD



Order No.

IMED-893

F-35 HEAD



Order No.

IMED-864

MODUS FOCUSED ELECTRODE



Order No.

IMET-890

F-60 HEAD



Order No.

IMET-891

F-80 HEAD



Order No.

IMED-886

REFILL AND DRAIN APPARATUS (FOR  
THE HEADS WITH COUPLING)



Order No.

IMED-886

REFILL AND DRAIN APPARATUS (FOR  
THE HEADS WITH NO COUPLING)



Order No.

IMED-892

F-L50 HEAD



Order No.

IMED-893


F-35 HEAD



## **> WHAT IS THE SHOCK WAVE THERAPY ?**


The administration of therapy using shock waves given outside of the body is a frequently preferred method in veterinary medicine, neurology, urology, cardiology, sports medicine, aesthetics, and especially in orthopedics and physical therapy. One of the advantages of this system is that it promotes faster tissue healing by increasing revascularization, collagen synthesis, and oxygenation in the application area. Furthermore, it develops less viscous tissue so that it is mechanically stronger.

With the non-invasive principle of the system, the desired treatment can be performed without the need for surgery.

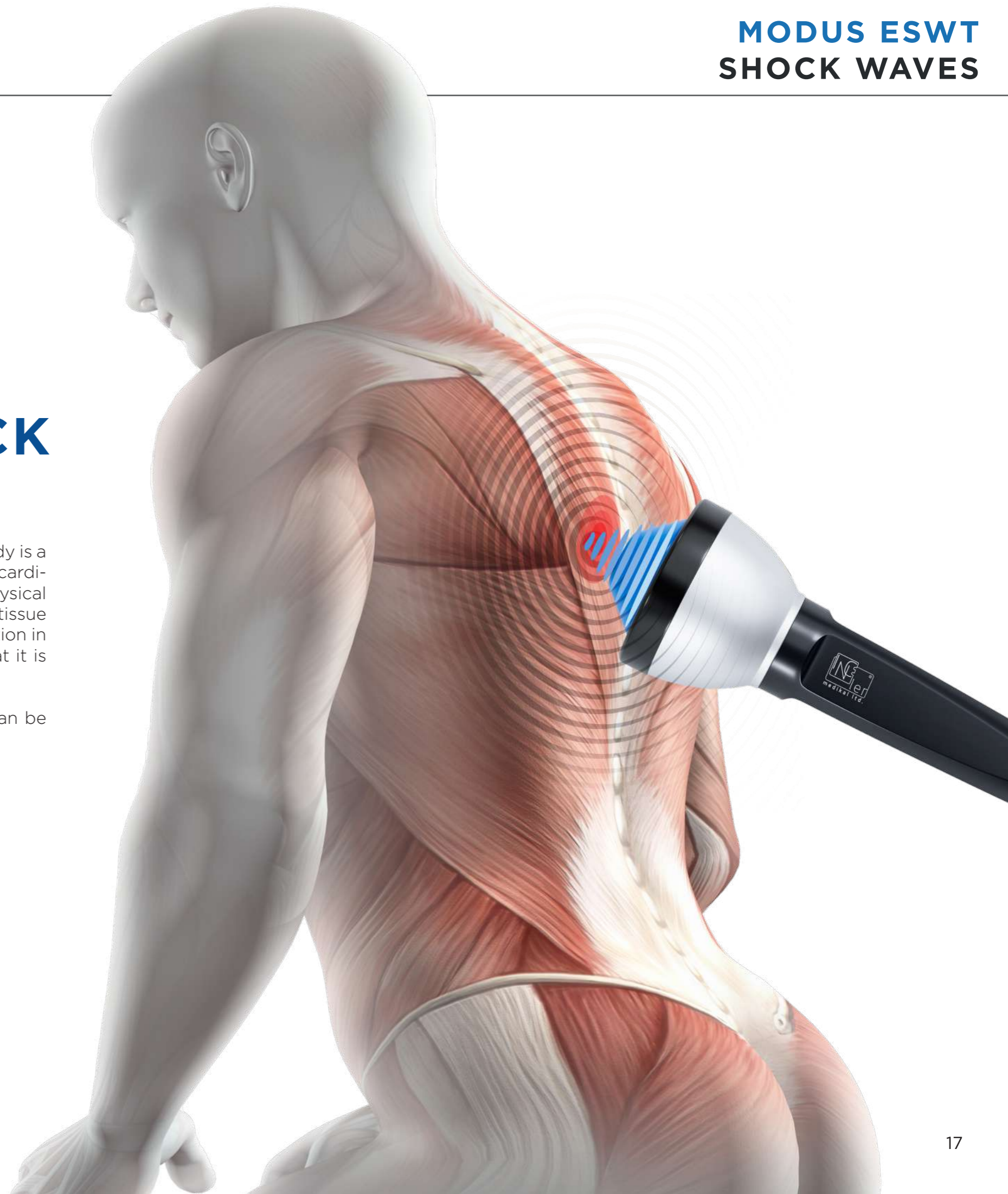



Orthopedics

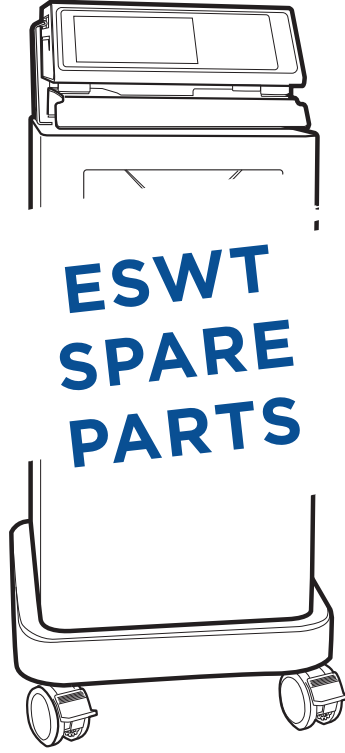

*inceler medikal*

**İNCELERMEDİKAL.COM**

**İnceler Medikal**  
İvedik OSB Mah. 1472. Cad.  
No: 120 Yenimahalle  
Ankara/Türkiye

tel: +90 549 255 3345  
fax: +90 312 255 3347  
sales@incelermedikal.com

